



Illustrative photo



ZEFIR

Instructions for use

Edition no.: 10/ 2025-11-12

Special needs stroller Zefir

This is a medical device.

Use it according to the instructions for use or label.

Purpose of the instruction

Instructions for use is designed for people who directly care of people with disabilities (assistants, parents, guardians), which will be transported in a stroller. It is also intended for doctors and physiotherapists operating the stroller.

ATTENTION! Before using the stroller, each User is obliged to read the instruction and act in accordance with its guidelines.

ATTENTION! In the event of any medical incident related to the Zefir stroller, the Manufacturer and the competent authority of the Member State in which the user or patient resides must be notified immediately, in writing, by registered letter or e-mail with acknowledgment of receipt.

Manufacturer's contact details:

- postal address : MyWam Sp. z o.o., Plebiscytowa 3, 41-100 Siemianowice Śląskie, Poland
- contact phone number : +48 32 733 11 31
- e-mail address : mywam@mywam.eu



Instructions for use version: **EN Edition 10/ 2025-11-12**

MyWam Sp. z o.o. reserves the right to make technical and commercial changes to the content of the Instructions for use without prior notice. Each change will be determined by the date the instructions were last updated. Soft copy of the operating instruction is available at www.mywam.eu

ATTENTION! Pictures given in the instructions for use are for reference only and can vary from the purchased model.

MyWam Sp. z o.o. has implemented a quality management system compliant with the ISO 13485 standard in the design, production and sale of class I medical devices. System certified by LL-C (Certification).



Table of contents

1. Purpose of Zefir stroller	5
2. Important safety instructions and rules.....	6
3. Plates.....	6
4. Symbols designation	7
5. General construction of the stroller.....	7
6. Configuration and use – stroller frame	8
6.1. Frame folding mechanism.....	8
6.2. Frame unfolding mechanism	8
6.3. Folding the stroller	9
6.4. Transport lock. Engaging and disengaging of the transport lock	9
6.5. Front swivel wheels. Front wheels direction lock, mounting, dismounting	9
6.6. Rear wheels. Mounting and dismounting of the rear wheel on the axle	10
6.7. Central brake. Engaging and disengaging of the central brake	11
6.8. Tilting foot.....	11
6.9. Push handle. Mechanism and angular adjustment.....	11
6.10. Footplate. Adjustment and use.....	12
6.11. Seat depth adjustment.....	12
6.12. Backrest height adjustment	13
6.13. Backrest recline adjustment	13
7. Standard equipment	14
7.1 . 5-point belt.....	14
7.2. Safety barrier.....	14
7.3. Abduction block (Wedge)	14
7.4. Abduction belt.....	15
8. Additional equipment	15
8.1. Shopping basket	15
8.2. Adjustable canopy with window.	15
8.3. Therapy table	16
8.4. Loin cushion	16
8.5. Headrest.....	16
8.6. Narrowing inserts.....	16
8.7. Travel bag	16
8.8. Raincover.....	17

8.9. Mosquito net	17
8.10. Sleeping bag	17
8.11. Sun umbrella	17
8.12. Winter gloves	17
9. Technical specification of the Zefir stroller	18
9.1. Material composition of the stroller	18
10. Delivery and tools	18
11. Assembly and disassembly of the stroller	19
12. Preparation for use	20
13. Using the stroller	20
13.1. Carrying the stroller	20
13.2. Getting on and off the stroller	20
13.3. Correct position in a stroller	21
13.4. Riding in a stroller	21
13.5. Moving on slopes	21
13.6. Overcoming thresholds or curbs	21
13.7. Transport in a car and other vehicles	21
13.8. Cleaning and care	21
14. Storage of the stroller	22
15. Service and periodic inspections	22
16. Disposal	22
17. Compliance with requirements for medical devices	23
18. List of authorized service points	24
19. Warranty card	25
20. Notes	26

1. Purpose of Zefir stroller

The special needs stroller Zefir is a medical device class I indicated for use by disabled person who, due to their health condition, cannot walk or sit down on their own. The product is intended for the transport and movement of the above persons in a sitting or lying position and for posture correction and this is its intended use.

Clinical indications include cerebral palsy, meningeal hernia, dystrophies, paresis, etc.

The product is intended for users 80-140cm height and up to 40 kg of weight. The product should be used in accordance with the recommendations of the attending physician or physiotherapist

Contraindications:

In case of contraindication to the sitting position and/or passive upright positioning, the adjustment of the stroller should be consulted with a physiotherapist or attending physician.

▪ Safety

2. Important safety instructions and rules

- It is forbidden to use the stroller in a way for which the stroller is not dedicated.
- The product contains small parts that can be dangerous if swallowed by a child!
- It is forbidden to put additional weight on the stroller, especially when it comes to placing an additional weight on the back of the backrest, e.g. a backpack or a bag. This may cause the Zefir stroller to tip over and injure the user! It does not apply to original accessories dedicated to the Zefir stroller that do not exceed the permissible load. All additional elements attached to the stroller affect its stability.
- ATTENTION!** While using, folding and adjusting of stroller mechanisms, there may be a risk of entrapment &/or squeezing some parts of the user body/caregiver in wholes/gaps between parts. All stroller adjustments must be carried out without any load on it (user load).

ATTENTION! Each time the stroller is used, the steps listed below must be performed and the stroller must not be used in the event of any irregularities. In such a situation, it is recommended to urgently contact the Seller and/or the Manufacturer for advice.

- ATTENTION!** When driving through doors or similar restrictions restricting the direction of travel of the stroller, make sure that there is adequate space on both sides of the stroller to prevent pinching of hands or other body parts and to prevent damage to the stroller.
- ATTENTION!** Temperature of the stroller may rise when sun exposed.

ATTENTION!



- Never leave user unattended in a stroller!
- This product is not suitable for running, skating or roller-skating
- This stroller is not intended as a seat in a motor vehicle
- Beware of gaps and moving parts – risk of jamming
- Do not carry the stroller with the user inside it.
- The stroller is packed in a plastic bag. Keep it away from infants and young children.


- ATTENTION!** It is forbidden to use the stroller if the technical condition of the stroller or a single element of the stroller is bad or the user has any doubts about it. In such a situation, it is recommended to urgently contact the Seller for advice.

3. Plates

Identification plate



Warning plate

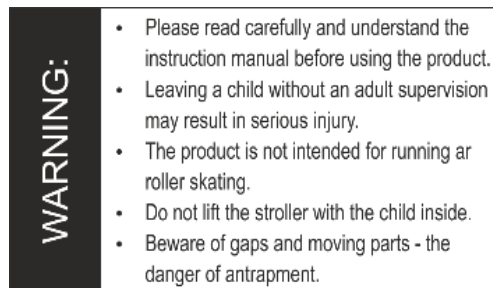

4. Symbols designation

	Product serial number		The manufacturer's name
	Permitted total user weight		The product is a medical device
	Date of manufacture		DESCRIPTION "WARNING!" and the symbol in this instruction are used to enhance the reader's attention to the content marked with this sign. Non-compliance with content marked with this description may endanger the health and safety of the user and accompanying persons, and may also damage the product.
	Conformity assessment mark according to applicable European regulations and directives		
	Warning about proceeding in accordance with the product Instructions for use		

Product

5. General construction of the stroller

Standard equipment



- 1- Handles for pushing the strollerSee point no.: 6.9
- 2- Backrest height adjustment mechanismSee point no.:6.12
- 3- Upholstery
- 4- Mechanism for adjusting the backrest angle See point no.: 6.13
- 5- Frame folding mechanism See point no.: 6.1
- 6- Transport lock See point no.: 6.4
- 7- Central brakeSee point no.: 6.7
- 8- Tilting footSee point no.: 6.8
- 9- Rear wheelsSee point no.: 6.6
- 10- Front swivel wheelsSee point no.: 6.5
- 11- Front wheels directional lockSee point no.: 6.5
- 12- Adjustable and removable footrestSee point no.: 6.10
- 13- Seat depth adjustment mechanism..... See point no.: 6.11
- 14- Abduction belt *See point no.: 7.4
- 15- Abduction block (Wedge) *See point no.: 7.3
- 16- Safety barrierSee point no.: 7.2
- 17- 5-point beltSee point no.: 7.1

*) options to choose with standard equipment

• Additional equipment

- 18- Shopping basket.....See point no.: 8.1
- 19- Travel bagSee point no.: 8.7
- 20- Sleeping bagSee point no.: 8.10
- 21- Narrowing insertsSee point no.: 8.6
- 22- HeadrestSee point no.: 8.5
- 23- Sun umbrella.....See point no.: 8.11
- 24- Adjustable canopy with window
.....See point no.: 8.2
- 25- RaincoverSee point no.: 8.8
- 26- Mosquito netSee point no.: 8.9
- 27- Loin cushionSee point no.: 8.4
- 28- Therapy tableSee point no.: 8.3
- 29- Winter glovesSee point no.: 8.12



6. Configuration and use – stroller frame

6.1. Frame folding mechanism

- left and right joint – element allowing to fold the upper part of the frame in order to minimize its size,
- frame folding mechanism with transport lock.



6.2. Frame unfolding mechanism

Place the frame flat on the surface and disengage the transport lock (Pict. 6.2.1). Next pull the backrest into the vertical position (Pict. 6.2.2).



Pict. 6.2.1



Pict. 6.2.2

In order to set up the stroller, carefully pull the handle up simultaneously setting the stroller in desired position till frame folding mechanism engages. Stroller frame will be correctly unfolded only once the frame folding mechanism

will be engaged (Pict. 6.2.3). This will be confirmed with a clear “click” sound.



Pict. 6.2.3



Pict. 6.2.4

6.3. Folding the stroller

It is important to make sure that the front wheels are unlocked and in the correct position (Fig. 6.3.2)

- **Folding:** unlock the frame folding mechanism by pulling up the frame folding locks simultaneously on the right and left side of the frame till its disengagement (Pict. 6.3.1). Fold the frame forwards, in direction of the front wheels mount. Engage the transport lock.

ATTENTION! In order to fold the frame, set the red clamp in low position.



Pict. 6.3.1



Pict. 6.3.2

6.4. Transport lock. Engaging and disengaging of the transport lock

- **Engaging the transport lock:** fold the frame, place the transport lock hole on the locking pin so that these two elements tightly overlap each other (Pict. 6.4).

- **Disengaging the transport lock:** before unfolding the frame, take the lock off the locking pin and then proceed with unfolding of the stroller frame.



Pict. 6.4

6.5. Front swivel wheels. Front wheels direction lock, mounting, dismounting

Two types of wheels can be used in the stroller: universal wheels (Fig. 6.5.1.) and slim wheels (Fig. 6.5.2).



Pict. 6.5.1

Universal wheels


Pict. 6.5.2

Slim wheels


Pict. 6.5.3

Quick release clamp

- **Engaging of the direction lock:** turn the lock knob right (in the direction of travel) - see Pict. 6.5.4.
- **Disengaging of the direction lock:** turn the lock knob opposite side – see Pict. 6.5.4.
- **Dismounting of the front wheels:** press the quick release clamp – see Pict. 6.5.3. and pull the wheel down.



Pict. 6.5.4

Direction lock knob

- **Mounting of the front wheels:** slide the wheel from below. Correct placement of the wheel will be confirmed with the clear “click” sound.

6.6. Rear wheels. Mounting and dismounting of the rear wheel on the axle

Two types of wheels can be used in the stroller: universal wheels (Pict. 6.6.1.) and slim wheels (Pict. 6.6.2).

- **Mounting the rear wheel on the axle:** press the quick release clamp located in the middle of the wheel (Pict. 6.6.3) and mount the wheel on the axle. Correct placement of the wheel will be confirmed with the clear “click” sound.



Pict. 6.6.1

- **Dismounting the rear wheel on the axle:** press the quick release clamp located in the middle of the wheel (Pict. 6.6.3) and slide the wheel out of its axle (Pict. 6.6.4).



Pict. 6.6.2



Pict. 6.6.3



Pict. 6.6.4

6.7. Central brake. Engaging and disengaging of the central brake

The central brake pedal (Pict. 6.7), which when pressed blocks the two rear wheels of the stroller directly, is located in the middle of the crossbeam, at the rear of the stroller frame!

- **Engaging:** press the brake pedal down towards the ground until it clicks. The correct activation of the brake is indicated by a clear "click" sound from the central brake.
- **Disengaging:** To disengage the central brake, pull the brake pedal up.



Pict. 6.7

6.8. Tilting foot

Using the tilting foot makes it easier for the stroller to climb onto small steps or curbs (Pict. 6.8).



Pict. 6.8

6.9. Push handle. Mechanism and angular adjustment.



Pict. 6.9.1

- **Adjustment mechanism.** These are two buttons on the outside of the frame, located on the joints of the handle to push the stroller. (Pict. 6.9.2).
- **Angular adjustment:** press both buttons of the angle

adjustment mechanism of the push handle and change its

angle. Once set in desired position, release the button to automatically lock the push handle.



Pict. 6.9.2

6.10. Footplate. Adjustment and use

- **Mounting:** slide the footrest perpendicularly onto the locking pins located at few different levels in the lower section of the frame. The clear „click” sound confirms correct placement of the footplate mounting hooks. Once mounted, make sure that the footplate it is mounted correctly. In order to check it, energetically pull up the footrest – if mounted properly it will remain in its position.



Pct. 6.10.1

Minimum position


Pct. 6.10.2

Maximum position


- **Dismounting:** simultaneously press the two mounting hooks located at both sides of the footplate and release the footplate. It is forbidden to use the stroller without the footrest.

- **Use:** It is recommended to position the footrest at such a height that the user's feet can rest freely on it with their entire surface, while the user's thighs are parallel to the seat surface (Pict. 6.10.1 i 6.10.2). There are straps in the footrest plate to secure the user's feet. It is recommended to put them on when using the stroller. Normally the footrest plate is in a horizontal position. However, during transport, storage, rehabilitation, getting on and off the Zefir stroller, the footrest plate must be dismounted.



- **Folding the footrest:** It is possible to fold the footrest plate.

Remark! During transport, storage, rehabilitation, getting in and out of the stroller, the footrest plate must be folded.

6.11. Seat depth adjustment

Thanks to the possibility of unfolding its front part, we can extend the seating area by 5cm.

Adjustment: press both buttons located on the sides of the front part of the seat simultaneously and bend the adjustable element to obtain the position of the main part of the seat (Pict. 6.11.1). To fold the seat, press the buttons located on the sides of the front part of the seat simultaneously and fold the adjustable element down or under the main part of the seat (Pict. 6.11.2).



Pict. 6.11.1

Seat depth adjustment - maximum


Pict. 6.11.2

Seat depth adjustment - minimum


6.12. Backrest height adjustment

Adjustment: detach the upper part of the backrest upholstery, insert the extensions (Pict. 6.12.2) into the holes of the backrest tubes and fasten the upper part of the upholstery again.



Pict. 6.12.1

Backrest height adjustment – max/min



Pict. 6.12.2

Extensions

6.13. Backrest recline adjustment



Pict. 6.13.1

*Backrest recline angle – sitting position
(basic)*



Pict. 6.13.2

*Backrest recline angle - semi-recumbent position
(reclined)*

• **Adjustment:** slide down the backrest mounts located on the sides of the frame and tilt the backrest from the basic position to the reclined position (Pict. 6.13.3). In order to lock the reclined position slide down the backrest mount and press it afterwards



Pict. 6.13.3

Lowering the backrest mount and reclining

7. Standard equipment

7.1 . 5-point belt

- **Adjustment:** move the adjusters on the straps up or down and move the protectors (Pict. 7.1.1 & 7.1.2). To adjust the height of the attachment of straps to the backrest, the harness attachment should be moved from the low attachment to the high attachment or conversely.
- **Use:** put the belt on the right and left user's arm and then attach both belts to the belt buckle that is fastened on the seat. To remove the five-point belt from user, press the button on the buckle and remove the belt from the user's left and right side.



Pict. 7.1.1



Pict. 7.1.2

5-point belt adjustment - max 5-point belt adjustment - min

ATTENTION! Every backrest angle adjustment requires readjusting of the length of the 5point belt.

7.2. Safety barrier



Pict. 7.2

Safety barrier

- **Mounting:** to be mounted on the stroller frame by using plastic latches, which look the same on both ends. After applying both latches at the appropriate height, press the them so that they completely cover the frame.

7.3. Abduction block (Wedge)

Abduction block should be placed between user's thighs.

- Mounting: attach it with Velcro on the front or back of the apron (Pict. 7.3.1 i Pict. 7.3.2),) and then wrap it with a band to strengthen it.



Pict. 7.3.1

Abduction block - max



Pict. 7.3.2

Abduction block - min

7.4. Abduction belt

Abduction belt is used for additional pelvis stabilization and prevents the user from slipping out of the stroller (Pict. 7.4). Wrap the user's thighs so that his pelvis adheres as close as possible to the backrest. This helps keep a healthier and more comfortable position of the user in the stroller.

- **Adjustment:** move the adjusters on the front belt straps up or down.
- **Use:** To put on abduction belt user should be seated on the straps, then wrap each user's leg and fasten both belts to the front of the two tapes. In order to remove belts, unfasten the belt buckles of the two front straps.



Pict. 7.4 Abduction belt

8. Additional equipment

8.1. Shopping basket

Shopping basket is located under the seat (Pict. 8.1). Its maximum load capacity is 2 kg.



Pict. 8.1

Shopping basket

8.2. Adjustable canopy with window.

- **Canopy elements:** frame (front and rear arch), canopy adjustment mechanism (right and left), two snaps with red buttons fastening the canopy to the frame, canopy upholstery, and a window with a mosquito

net with a zipper.

- **Mounting:** slide down both canopy latches on the frame mounts till clear "click" sound.

- **Dismounting:** press the two red buttons on the canopy fastening latches (located on the inside), then slide the canopy fasteners out of the mounts on the frame.

Use: In order to change the size of the canopy window area, hold the front arch and the rear arch of the canopy frame and bring the arches closer or closer to each other to set the desired position. The mechanism adjusting the size of the canopy surface locks automatically. The canopy has a removable window with a mosquito net, ensuring free air circulation. (Pict. 8.2.1). The window opens and closes with a zipper.



Pict. 8.2.1

Canopy – max. setting



Pict. 8.2.2

Canopy – min. setting

8.3. Therapy table

Maximum load capacity of the table is 3kg (Pict. 8.3).

- **Mounting:** slide the open fastening clamps onto the main element of the frame, and then tighten the clamps until they are completely closed.
- **Dismounting:** unfasten the clamps and remove the table from the stroller frame.



Pict. 8.3
Therapy table

8.4. Loin cushion

Loin cushion increases the comfort of using the stroller by softening the backrest in the lumbar section (Pict. 8.4).

- **Mounting:** attach it to the backrest with velcro.

8.5. Headrest

Headrest performs both softening and stabilizing functions (Pict. 8.5).

- **Headrest height adjustment:** attach the it with velcro at the correct height of the backrest. Then, for better stability, drag the headrest straps fitted with velcro through the backrest openings and attach them behind the backrest.

8.6. Narrowing inserts

Narrowing inserts are used to adjust the inside dimensions of the seat width (Pict. 8.6).

- **Mounting:** attaching them with the longest vertical edge of the pad to the outer edges of the backrest and to the side elements of the stroller frame using a velcro strap.



Pict. 8.4 Loin cushion



Pict. 8.5 Headrest



Pict. 8.6 Narrowing inserts

8.7. Travel bag

Elegant travel bag is mounted on two handles, located on the inside of the handle to push the stroller (Pict. 8.7). Maximum load capacity of the bag is 2kg.

ATTENTION! Exceeding the maximum bag load has a negative effect on the stability of the stroller and there is a danger for the user and other people!



Pict. 8.7
Travel bag

8.8. Raincover



Pict. 8.8
Raincover

This accessory requires to be installed to other optional equipment - an adjustable canopy with window.

- **Mounting:** put it on top of the entire canopy, and then cover the seat and footrest (Pict. 8.8).

8.9. Mosquito net

This accessory requires to be installed to other optional equipment - an adjustable canopy with window (Pict. 8.9).

- **Mounting:** put it on top of the entire canopy, and then cover the seat and footrest. When mounted properly, attach it to the frame and footrest.



Pict. 8.9 Mosquito net

8.10. Sleeping bag



Pict. 8.10 Sleeping bag

This accessory requires to be installed to other optional equipment – safety barrier (Pict. 8.10).

- **Mounting:** put it on the footrest from below, and attach it from above through the safety barrier. Then fasten it with snaps to the frame of the stroller.

8.11. Sun umbrella

Mounting: loosen the clamp by unscrewing the thread located on the side, then slide it on, ensuring that the rubber washer remains in place. After choosing the desired height, tighten the thread until the clamp on the stroller frame is fixed properly (Pict. 8.11).



Pict. 8.11 Sun umbrella

8.12. Winter gloves

Place on the push handle with the help of the sewn-in magnets (see Pict. 8.12).



Pict. 8.12 Winter gloves

9. Technical specification of the Zefir stroller

	Universal wheels	Slim wheels
length of the stroller	95 – 113 cm	91-111 cm
the width of the stroller	62 cm	60 cm
the height of the stroller	102 – 122cm	100-118 cm
total weight	13 – 17 kg	13-16 kg
effective seat depth	24 / 30 cm	
seat width	29* / 38 cm	
backrest height	50 / 60 cm	
backrest width	29* / 37 cm	
distance between footrest and the seat / number of adjustments	20, 30, 40 cm / 3 positions	
length of the stroller, folded up for transport	94 cm	90 cm
width of the stroller, folded up for transport	62 cm	60 cm
height of the stroller, folded up for transport	35 cm	35 cm
capacity of the shopping basket	2 kg	
carrying capacity of a travel bag	2 kg	
Maximum user's height	40 kg	
Maximum user's weight	80-140 cm	
seat plane angle / number of adjustments	20° / 1 position	
backrest angle / number of adjustments	25° i 57° / 2 positions	
footrest angle / number of adjustment	115° / 1 position	

COMMENTS: *) values shown with inserts

Measurement value tolerance +/- 20mm, +/- 1.5kg, +/- 1 °

The manufacturer reserves the right to introduce changes to the technical parameters.

9.1. Material composition of the stroller

Powder coated steel and aluminium, fabrics 100% polyester and plastic (PP, PCV, ABS, PA, PE, PU).

10. Delivery and tools

The original packaging should contain the following items:

- Folded stroller frame with integrated seat, backrest and wheels, footrest, safety barrier and wedge or abduction and stabilizing straps
- Additional equipment (according to the order)
- Instructions for use with warranty card.

Before use:

- a. Make sure that the package has not been damaged. If it is found that the package is damaged, then at the presence of the courier, unpack the goods and check whether none of the stroller elements hasn't been damaged.
- b. On the basis of the specification check the package content (in case of any missing elements contact your Seller)

No tools are required to assemble and disassemble the Zefir stroller.

11. Assembly and disassembly of the stroller

- **Assembling the stroller.**

The product pre-packed in the original box is already pre-assembled. Complete the installation by completing the following steps:

- unfold the stroller frame by adjusting the parameters
- engage / disengage the front direction lock
- adjust the height and angle of the handle to push the stroller
- set the angle of the backrest
- adjust the length of the 5-point belts,
- put on a footrest at the desired height,
- put on a safety barrier or a therapy table,
- mount appropriate equipment and adjust parameters, for example:
 - put narrowing inserts,
 - remove / install the extension on the backrest
 - set the seat depth,
 - adjust the abduction belt
 - adjust the canopy with the window
 - install and adjust the abduction block (Wedge)
 - put on the travel bag,
 - put on the shopping basket
 - sleeping bag,
 - mosquito net,
 - rain cover,
 - sun umbrella
 - winter gloves.

- **Disassembling the stroller**

In order to disassemble the stroller:

- Dismount all the additional equipment,
- Dismount the footrest,

- Dismount the safety barrier / therapy table
- Lock the front wheels,
- Engage the central brake on the rear wheels
- Fold the stroller closing it with the transport lock

12. Preparation for use

Before using the stroller:

- Check that the stroller frame has been folded out correctly and that the folding mechanism has been properly closed (locked).
- Correctly make all regulations and adjust the stroller to the user's needs, including setting the appropriate depth of the therapeutic seat and the height of the backrest, appropriate backrest recline, set the appropriate height of the footrest, adjust the appropriate tension of the 5-point belts. Check if all adjustments are properly seated and secured.
- Check if the stroller upholstery is dry.
- Secure the user with all available stabilizing elements, such as: 5-point belts, narrowing inserts, abduction belt, abduction block (wedge), footrest, footrest tapes, safety barrier or therapy table. This is to prevent the user from falling out and sliding out of the stroller and to prevent the user's feet from falling off the footrest. Make sure that all adjustments made earlier are suitable and safe for the user.
- Make sure that the user is not placed in the stroller without clothing that protects the body at the points in contact with the stroller's upholstery. This is to prevent chafing of the user's body and to keep the user's hygiene.
- In the case of users who have problems with incontinence, a special oilcloth should be placed on the stroller seat, which can be purchased in pharmacies or medical stores.
- In the case of using any additional aids, such as vents, drains, bags, catheters, etc., make sure that they are well protected against falling out, pulling out, tightening and not interfering with the use of the stroller.

13. Using the stroller

13.1. Carrying the stroller

Appropriate way to carry a stroller is to drive it using its wheels. In cases where this is not possible (e.g. to transport the stroller in a car), for safety reasons, the stroller can be lifted and moved with both hands after dismantling all additional equipment to reduce the carrying weight and to fold the stroller properly.

Be careful when moving or lifting the stroller's frame. Frame should be carried once it is folded and secured with the engaged transport lock.

Moving of the stroller with the help of any of its accessories is forbidden – possibility of its damage.

13.2. Getting on and off the stroller

During getting on and off the stroller, place it as close as possible to a stable and safe place where the user will change. After making sure that the rear brakes are applied and the footrest has been removed to avoid standing on it, the user can be placed in the stroller.

Do not stand on the footplate, as there is a risk of tipping over the stroller, injuring the user and damaging the product. Make sure that when getting in or out of the stroller, the user and any other person should not stand on the footrest, as this may damage it and injure the user!

13.3. Correct position in a stroller

For safety and health reasons, it is recommended to place user's pelvis and back as close as possible to the backrest to support the spine properly (optimal support of all parts of the spine so that there is no free space).

Thighs should be supported on the largest surface, along their entire length, which is used to adjust the depth of the seat of the stroller and the feet should rest on the footrest plate, which prevents slipping, maintains optimal tension of the leg muscles and carries part of the weight from higher parts of the body.

13.4. Riding in a stroller

When using the stroller, always hold it by the handles to push to prevent the stroller from rolling. Release the rear brakes while holding the push handle and push the stroller in the desired direction.

13.5. Moving on slopes

Always position the stroller forwards when climbing or driveways. Limit the speed on sloping surfaces and move as slowly as possible. The stroller assistant (parent, guardian) must have sufficient strength to control the stability of the stroller.

Always remember that the stroller's centre of gravity changes when it overcomes unevenness or slopes on the ground. In this case, be extremely careful and secure the stroller. When using ramps, make sure that the ramp is suitable, that the wheels of the stroller are high enough to overcome the ramp, that the frame of the stroller will not touch the ramp, and that a second person should be provided as there is a risk of the stroller tipping over, damaging the product and injuring the user.

13.6. Overcoming thresholds or curbs

When overcoming small obstacles (small thresholds, small curbs, etc.), block the tilting foot with the foot and the handle to push the stroller downwards and towards you to raise the front wheels of the stroller and carefully drive the front wheels over the obstacle. Then slowly pull the push handle upwards away from the ground to raise the rear wheels of the stroller and carefully drive the rear wheels onto the obstacle.

13.7. Transport in a car and other vehicles

- Remove all accessories from the stroller and store them in a safe place. Then dismantle the stroller and place it in the trunk. If the dismantled stroller is not in the trunk of the car separated from the passenger compartment, then the stroller frame must be securely attached to the vehicle. For this purpose, unused seat belts may be used.
- In order to transport a stroller in other vehicles, such as trains, trams, buses, planes, ships, etc., please ask the persons responsible for transport in these means of transport about the conditions of safe travel to the destination with a stroller and the appropriate requirements.



The Zefir stroller should not be used as a seat to transport an user in a vehicle.

13.8. Cleaning and care

The upholstery can be washed in a washing machine at temperatures up to 40°C.

The upholstery should be dried on the stroller structure. All stains should be cleaned immediately after.

ATTENTION! Do not clean the upholstery with the use of any chemical compounds composed of solvents or bleaches. Manufacturer does not take any responsibility for damages caused by the use of inappropriate cleaning products. Make sure the product is always dry before use.

When using detergents designated for this type of material follow its instruction manual.



Wash at 40°



Do not bleach



Do not tumble dry



Do not iron

14. Storage of the stroller

The stroller should be stored in a clean, dry and well ventilated room, preferably at 0-35 Celsius degrees. In case of a long-term storage it is necessary to follow each steps described in point no.: 15 on a monthly basis.

15. Service and periodic inspections

For safe and trouble-free use of the Zefir stroller, it must be regularly and thoroughly inspected and serviced in designated service points.

In a situation where the stroller user is not able to carry out the activities himself, he should return the stroller to a specialized stroller service point or directly to the Manufacturer, who will properly review and will perform necessary service activities. These activities may be charged extra.

All repairs should be carried out in specialized, authorized Services (see point 18) or directly at the Manufacturer. Only original spare parts that are supplied or recommended by the Manufacturer may be used in the stroller.

Service should be performed:

- **Each time before using the stroller**

- Check the efficiency of the central brake - After activating the central brake, both rear wheels of the stroller should be locked and cannot rotate.
- Check the attachment of structural elements - It is necessary to check all mechanisms of the stroller and its movable elements, in particular, pay attention to: correct mounting of the seat with the backrest to the frame of the stroller, proper mounting of all the 4 wheels, proper work of the 5-point belt, correct functioning of the seat angle adjustment mechanism, correct functioning of the backrest recline mechanism and proper mounting of the footrest.
- Check the condition of the tires - Check that the tires on all four wheels are not damaged or worn.
- Inspection of solid frame parts, seat with backrest and footrest frame - Visually inspect the condition of the structural profiles. These elements must not be damaged or cracked.

- **Every one month**

- Cleaning of the stroller
- Visual check of all the moving elements, bearings, adjusting mechanisms – lubricate with silicone grease.
- Check all bolted connections for looseness.

16. Disposal

Disposal and recycling: in accordance with the regulations in force in the given country.

17. Compliance with requirements for medical devices

Class I medical device in accordance with the Regulation of the European Parliament and of the Council (EU) 2017/745 of April 5, 2017 on medical devices. The manufacturer issued an appropriate declaration of conformity for the product, which is available at www.mywam.eu

18. List of authorized service points

Poland
MyWam Sp. z o.o.
Wiosny Ludów 91, 40-373 Katowice, Poland
phone: +48 32 733 11 31
Website: www.mywam.eu
Email address: mywam@mywam.eu

19. Warranty card

Warranty Card Series and number:

Date of purchase: - -

Seller stamp and signature:

.....
.....
.....

Warranty conditions

- The Manufacturer (Guarantor) provides the Buyer with a warranty for the purchased stroller for a period of 24 months, counting from the date of purchase.
- The buyer complaining about the stroller should deliver the stroller to the Seller in a clean package (free from stains, additional equipment not belonging to the stroller, moisture, hair, etc.) along with the stroller elements protected against damage.
- Damage caused by insufficient protection of the forklift shipment and its components is not subject to warranty repairs.
- The complainant (Buyer) may not demand that the Manufacturer repair the damage in the event of damage to the shipment if he incorrectly labels or incorrectly insures the shipment of the forklift and its components.
- The warranty period in the case of accepting the complaint is extended by the time of repair of the stroller, which is counted from the moment of handing over the stroller to the Seller under warranty to the day when the stroller was ready for collection.
- The activities provided for in the instruction, which the stroller user is obliged to perform on his own and at his own expense, cannot be performed as part of the warranty repair.
- Failure to use the stroller in accordance with the intended use, inconsistent with the instructions for use or inconsistent with the principles of proper operation will void the warranty.
- Any structural changes to the truck void the warranty (in the area of change or in its entirety).
- Damage to the forklift caused by: improper operation, improper maintenance, improper storage of the stroller, mechanical damage and natural abrasion and wear of parts (e.g. upholstery, wheels) are not covered by the warranty.
- The manufacturer warns against warranty repairs in unauthorized service facilities. Finding such a repair or unauthorized modification of the stroller will void the warranty.
- The warranty card is valid provided that it has: the date of sale confirmed by the Seller's company stamp, the Seller's signature, the type and serial number of the truck.
- The warranty card is valid provided that it does not contain deletions or corrections.
- The manufacturer of the forklift conducts a paid post-warranty service at its headquarters.
- The manufacturer specifies the expected service life of the product at 4 years. Further use after this period is possible, but taking into account the care for safety and comfort of use, we recommend that the manufacturer or authorized service center carry out a technical inspection. The inspection service is payable and includes an assessment of the technical condition of the device, safety verification and any necessary maintenance activities.
- In matters not covered by this warranty card, the provisions of the Civil Code (Poland) shall apply.
- All matters related to claims under the warranty or warranty of the stroller should be directed to the Seller.

20. Notes